



i

Flexible Sinuslift Küretten

Optimales Feingefühl beim lateralen Sinuslift.


Die versehentliche Perforation der Sinusmembran ist eine Komplikation, die häufig während der Anhebung des Sinusbodens (Sinuslift, SBE) auftritt. Das neue SCHWERT Sinuslift-Instrumentarium erhöht durch stumpfe, flexible Arbeitsenden mit „Swing-Effekt“ die Taktilität und das operative Feingefühl bei der direkten Sinusbodenelevation. Durch die spezielle Formgebung und Beschaffenheit des Edelstahls zwischen Griff und aktivem Part wird leichter Druck abgefedert und dem Operateur eine direkte Rückmeldung vermittelt. Gleichzeitig wird die Perforationsgefahr der Schneider'schen Membran auf ein Minimum reduziert. Bereits bei der Entwicklung wurde auf eine besonders praktikable Handhabung geachtet. Durch die verschieden abgewinkelten Enden wird die Navigation auch in schwer zugänglichen Bereichen erleichtert. Die Instrumentenlinie umfasst sechs Kürettenfiguren um der jeweiligen OP-Situation optimal gerecht zu werden. Ergänzt werden diese durch Stopfer-/ und Löffelinstrument zur Einbringung des Augmentats.

Flexible sinus lift curettes

An optimum sensitivity during the laterally lift procedure.

The perforation of the Schneiderian membrane is a complication that often occurs during a maxillary sinus floor augmentation procedure. The new SCHWERT Sinus-Lift Instruments increase the sensitivity by blunt and flexible working tips with "Swing-Effect". Because of their special shapes and the characteristic of the stainless steel between the handle and the active part, light pressure is cushioned and gives the operator a direct and sensitive response. The risk to perforate the Schneiderian membrane is reduced to a minimum.

Already during their development the focus was a practicable handling. Because of varying angulations and shapes of the active parts they are easy to navigate also in difficult areas. The line of curettes consists of six different figures for every situation. Spoon and condenser complement the set to place the bone graft.

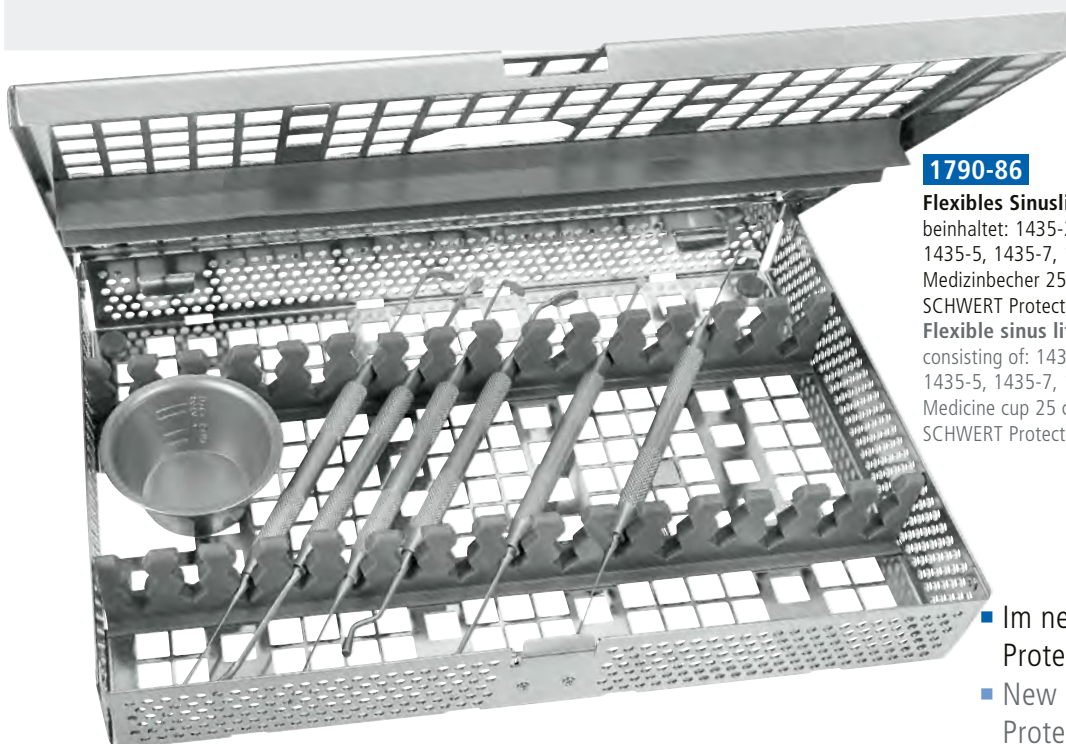

1790-86

Flexibles Sinuslift Küretten-Set, beinhaltet: 1435-2, 1435-3, 1435-4, 1435-5, 1435-7, 1435-8, 2616-25 Medizinbecher 25 ccm, SCHWERT Protector-Korb 2603-01
Flexible sinus lift curettes set, consisting of: 1435-2, 1435-3, 1435-4, 1435-5, 1435-7, 1435-8, 2616-25 Medicine cup 25 ccm, SCHWERT Protector-Basket 2603-01

- Im neuen SCHWERT Protector-Korb
- New SCHWERT Protector-Basket



SCHWERT
SWORD • ESPADA
EPEE • SPADA



1435-1

FIG. 1
20,4 cm



1435-2

FIG. 2
20,7 cm



1435-3

FIG. 3
20,0 cm



1435-4

FIG. 4
19,5 cm



1435-7

FIG. 5
19,0 cm



1435-8

FIG. 8
20,0 cm



1435-5

17,6 cm
Löffelinstrument
Spoon instrument